
beurer

IPL 7800 Cool Pro

- | | |
|---|--|
| DE IPL Haarentfernungsgerät
Gebrauchsanweisung4 | PL Depilator IPL
Instrukcja obsługi 143 |
| EN IPL hair removal device
Instructions for use.....24 | NL IPL-ontharingsapparaat
Gebruiksaanwijzing 163 |
| FR Epilateur à lumière pulsée IPL
Mode d'emploi.....43 | DA IPL hårfjerningsapparat
Betjeningsvejledning..... 183 |
| ES Depiladora de luz pulsada
Manual de instrucciones63 | SV IPL hårborttagningsapparat
Bruksanvisning.....202 |
| IT Epilatore IPL
Istruzioni per l'uso83 | NO IPL hårfjerningsapparat
Bruksanvisning.....221 |
| TR IPL epilasyon cihazı
Kullanım kılavuzu..... 103 | FI IPL ihokarvojen poistolaitte
Käyttöohje.....241 |
| RU Эпилятор IPL
Инструкция по применению.. 122 | |

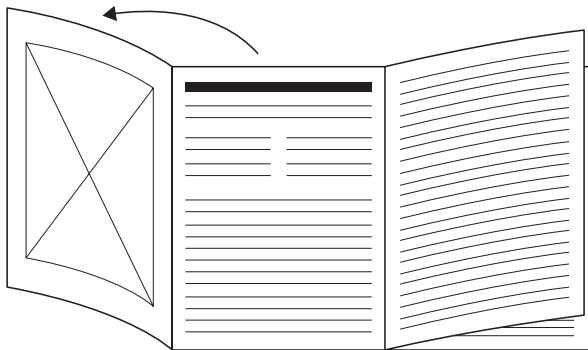
- DE** Klappen Sie vor dem Lesen der Gebrauchsanweisung die Seite 3 aus.
- EN** Unfold page 3 before reading the instructions for use.
- FR** Dépliez la page 3 avant de lire le mode d'emploi.
- ES** Despliegue la página 3 antes de leer las instrucciones de uso.
- IT** Prima di leggere le istruzioni per l'uso aprire la pagina 3.
- TR** Kullanım kılavuzunu okumadan önce 3. sayfayı açın.
- RU** Перед чтением инструкции по применению разложите страницу 3.
- PL** Przed przeczytaniem instrukcji obsługi otworzyć stronę 3.
- NL** Vouw pagina 3 uit om de gebruiksaanwijzing te kunnen lezen.
- DA** Fold side 3 ud, før du læser betjeningsvejledningen.
- SV** Vik ut sid. 3 innan du läser bruksanvisningen.
- NO** Åpne side 3 før du leser bruksanvisningen.
- FI** Käännä sivu 3 auki ennen käyttöohjeen lukemista.

1


2

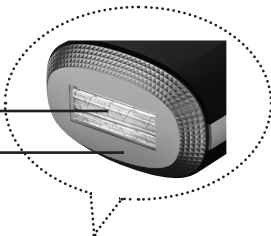3

4


5

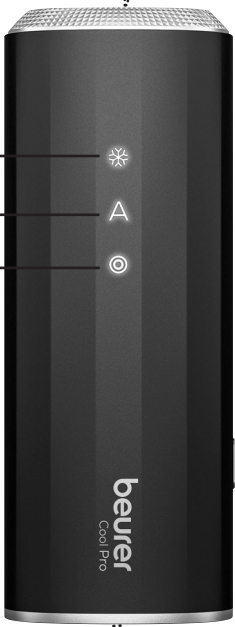7

8

9

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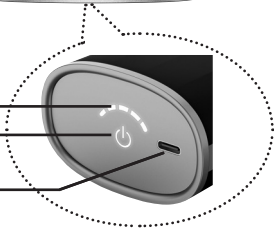Read these instructions for use carefully. Follow the warnings and safety notes. Keep these instructions for use for future reference. Make the instructions for use accessible to other users. If the device is passed on, provide the instructions for use to the next user as well.

TABLE OF CONTENTS

1. Included in delivery	26
2. Signs and symbols	27
3. Intended use	29
4. Warnings and safety notes	30
5. Device description	34
6. First steps	34
6.1 Preparing for use	34
6.2 Skin sensitivity test.....	35
7. Usage	36
7.1 Application methods	37
7.2 Special instructions for facial treatments	38
7.3 Application cycle	39
7.4 After treatment	40
8. Cleaning and maintenance	40
9. What if there are problems?	40
10. Disposal	41
11. Technical specifications	41
12. Warranty	42

⚠ WARNING

- The device and adapter are only intended for use in the home/private use and are not intended for commercial use.
- The device is only suitable for adults over the age of 18.
- Children must not play with the device.
- The device and adapter must not get wet.
- Check that the device and the mains cable are in perfect working order before each usage. In the event of damage, the device must not be used.
- Never look directly into the light emitted from the IPL lamp.
- Do not insert any objects into the device.
- Make sure that the device's vents are not covered.
- Do not use the device near to damp environments, such as bathtubs, wash basins, showers or other containers with water or other liquids – risk of electric shock!
- In order to avoid an electric shock, do not touch any internal components if the device is damaged.
- If the mains connection cable of this device is damaged, it must be disposed of. If it cannot be removed, the device must be disposed of.
- Never operate the device unsupervised. Supervise children around the device to ensure they do not play with it.

- Under no circumstances should you open or repair the device yourself, otherwise proper functionality can no longer be guaranteed. Failure to comply will invalidate the warranty.
- This device may be used by persons with limited physical, sensory or mental capabilities, provided they are supervised by a responsible person or they receive detailed instructions on how to use the device.

1. INCLUDED IN DELIVERY

Check that the exterior of the cardboard delivery packaging is intact and make sure that all contents are present. Before use, ensure that there is no visible damage to the device or accessories and that all packaging material has been removed. If you have any doubts, do not use the device and contact your retailer or the specified Customer Services address.

- 1 IPL hair removal device
- 1 mains adapter
- 1 these instructions for use

2. SIGNS AND SYMBOLS

The following symbols are used on the device, in these instructions for use, on the packaging and on the type plate for the device:

DANGER

Indicates an imminent danger. If it is not avoided, it will result in death or serious injury.

WARNING

Indicates a potentially impending danger. If it is not avoided, there is a risk of death or serious injury.

CAUTION

Indicates a potentially impending danger. If it is not avoided, slight or minor injuries may result.

NOTICE

Indicates a potentially harmful situation. If it is not avoided, the system or something in its vicinity may be damaged.

	Danger This device must not be used near to water or in water itself (e.g. in a wash basin, or in the shower or bathtub) – risk of electric shock!		CE labelling This product satisfies the requirements of the applicable European and national directives.
	Warning; Optical radiation		Polarity of DC power connection
	Keep away from eyes!		Do not use on tanned skin or after exposure to the sun
	Product information Note on important information		Read the instructions

	Disposal in accordance with the Waste Electrical and Electronic Equipment EC Directive (WEEE)		Separate the packaging elements and dispose of them in accordance with local regulations.
	Manufacturer		For indoor use only
	The products demonstrably meet the requirements of the Technical Regulations of the EAEU		Separate the product and packaging elements and dispose of them in accordance with local regulations.
	Marking to identify the packaging material. A = material abbreviation, B = material number: 1-7 = plastics, 20-22 = paper and cardboard		Protection class II device The device is double-insulated and, therefore, meets protection class 2
	Conformity assessed marking for the United Kingdom	ta=35°C	Nominal ambient temperature Maximum operating temperature
	Alternating current		Direct current The device is suitable for use with direct current only
	Protected against solid foreign objects, 12.5 mm in diameter and larger		Do not use the mains adapter if the plug is damaged.
	Switching mains adapter		Safety isolating transformer, short-circuit proof
	Energy efficiency level 6		Importer symbol

3. INTENDED USE

The device is intended for the removal of unwanted body hair in both women and men. The device must only be used on the **legs, armpits, bikini line, chest, stomach, back and face (only women may use it on the face, below the cheek line)**.

The device may only be used for the purpose for which it is designed and in the manner specified in the instructions for use. Any form of improper use can be dangerous. The manufacturer is not liable for damage resulting from improper or incorrect use.

IPL technology – how it works

The IPL 7800 Cool Pro is a light-based device for long-lasting hair removal for private use. The light-based hair removal process is a well-known and proven method. In over 15 years of clinical applications around the world, it has proven to be a safe and effective method for achieving long-lasting hair reduction. The additional cooling function can make usage more pleasant.

How does the light remove the hair?

Light-based hair removal is founded on the principle of selective photothermolysis, which uses optical energy to deactivate hair growth. To achieve such a thermal effect, the hair root must selectively absorb the light energy and convert it into heat.

This selectivity is achieved when the optical energy that is released into the tissue is mainly absorbed by the melanin pigment in the hair. The high-energy light is converted into thermal energy (heat), which ultimately deactivates hair growth.

What effect does the hair growth cycle have on light-based hair removal?

Anagen – growth phase

Catagen – transitional phase


Telogen – resting phase

Every hair on our body goes through the three phases of the hair growth cycle: Anagen, catagen and telogen. Only hairs that are in the anagen phase, or growth phase, respond to treatments with IPL.

Although the duration of a complete hair growth cycle differs from person to person and depends on the area of the body the hair is growing on, the cycle generally lasts 18–24 months.

Please note that at least one complete hair growth cycle is required in order to achieve long-lasting hair removal using the device.

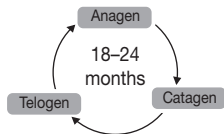4. WARNINGS AND SAFETY NOTES

⚠ DANGER

Keep the device and adapter away from water! The IPL 7800 Cool Pro is an electrical device. Do not place it in water or allow it to fall into water or any other liquid.

This can lead to severe electric shock!

- Do not use the device when taking a bath.
- Do not use the device when it is damp or wet.
- Do not reach for the device if it has fallen into water.
- Immediately unplug the mains adapter if it has fallen into water.
- Never touch the mains adapter or the device with wet hands!

⚠ WARNING

- Do not use the device when pregnant or breastfeeding.
- If you discover extreme reddening, blistering or burns on the skin, stop using the device immediately!
- Never attempt to release a light pulse into the air!
- On the basis of the “Hair and skin colour chart” on the back page of these instructions for use, determine which light intensities are best suited to the colour of your skin and your body hair and whether this method is suitable for you to use at all (if not, this is marked by an “x” in the chart).
- Do not use the device on naturally dark skin.
- Do not use the device with a contact gel – this may damage the device.
- Eye injuries (e.g. possible loss of vision) or skin injuries are possible if the instructions are not followed. Read and follow the instructions.
- If unexpected side effects occur within 48 hours of first use, e.g. itching, seek medical advice before further use.

The device works best on dark hair types or on hair with a high melanin content. Black and dark brown hair shows the best response and, although brown and light brown hair will also respond to the treatment, it usually requires more hair removal sessions. Red hair could produce a reaction. White, grey or blond hair does not usually respond to the device, although some users have noticed results after multiple hair removal sessions.

⚠ WARNING

Do not use the device under any circumstances if you suffer from one of the following ailments:

- If you suffer from skin cancer or have done in the past, or if you have had potential malignancies on the skin areas to be treated.
- If you suffer from epilepsy with flash sensitivity.
- If you suffer from collagen disorders, including the formation of keloids or poor tissue repair.
- If you suffer from herpes or psoriasis in the treatment area, unless you have consulted a doctor and received preventative treatment.
- If you are photosensitive or suffer from another form of light sensitivity, such as porphyria, polymorphic photodermatosis, solar urticaria, lupus, etc.
- If you suffer from a metabolic disorder, such as diabetes.
- If you suffer from a blood clotting disorder.

Do not use the device under any circumstances if you are taking any of the following medication:

- If you are taking medication that makes the skin more sensitive to light, including non-steroidal anti-inflammatory drugs (for example, aspirin, ibuprofen, paracetamol), tetracyclines, phenothiazines, thiazide diuretics, sulphonylurea, sulfonamide, DTIC, fluorouracil, vinblastine, griseofulvin, alpha hydroxy acids (AHAs), beta hydroxy acids (BHAs), Retin-A®, Accutane® and/or topical retinoids.
- If you are undergoing or recently underwent treatment with azelaic acid.
- If you have been on a course of steroids in the past 3 months.
- If you have had a skin exfoliation treatment or other skin smoothing procedure in the last 6–8 weeks.
- If you have undergone radiation or chemotherapy in the past 3 months.

Do not use the device under any circumstances in the following areas:

- On mucous membranes in the ear or nose area.
- On or next to an active implant, such as a pacemaker or a device for incontinence, insulin pump, etc.
- On the face above the cheek line, around the eyes, eyebrows or eyelashes, as this could result in severe damage to the eyes. 
- On nipples, areolae, inner labia, vagina, anus, nose and ears.
- Men must never use the device on the scrotum or face.
- On piercings or other metal objects, such as earrings or jewellery.
- On tattoos or permanent make-up in the area to be treated.
- On dark brown or black patches, such as large freckles, birthmarks, moles, blisters or warts.
- On eczema, psoriasis, lesions, open wounds or acute infections. Wait until the affected area has healed before using the device.
- On damaged or cracked skin or in areas at risk of malignancy.
- On an area of your body where you may like to have hair in the future.

- Do not use the device on an area of the body into which fillers or toxins have been injected (e.g. Botox).

This list is not exhaustive.

If you take prescription medicine, are monitored regularly due to an illness and/or are unsure whether the device is safe for you to use, consult your doctor or dermatologist before using the device.

Possible side effects:

If the device is used in accordance with these instructions for use, side effects and complications rarely occur in relation to the use of the device.

- After usage, it is possible that the skin will be slightly red, itchy or warm to the touch. These reactions are harmless and will quickly wear off.
- For a short time, you may experience a skin reaction similar to sunburn on the treated area. If this reaction does not wear off within three days, consult a doctor.
- Dry or irritated skin may occur as a result of shaving or a combination of shaving and light therapy. This reaction is harmless and will wear off within a few days. If you wish, you can cool the affected area with some ice or a damp cloth.
- If the skin continues to be too dry, you can apply an unscented moisturising cream 24 hours after using the device.
- In rare cases, the treated skin can become very red and swollen. This is usually because the intensity level you selected was too high and is more likely to occur in sensitive parts of the body. Reddening and swelling should subside within 2 to 7 days and should be treated by frequent application of ice. Gentle cleaning is OK, but you should avoid exposure to the sun.
- Discolouration of the skin is very rare and appears as darker or lighter patches on the affected skin area. It can be caused by using an intensity level that is too high for your skin tone. If the discolouration does not subside within two weeks, consult a doctor. Do not treat the discoloured areas again until the discolouration has fully disappeared.
- In very rare cases, burns or lesions to the skin can occur after using the device. The burn or lesion may take a few weeks to heal and, in extremely rare cases, you may be left with a visible scar.
- In very rare cases, skin infections and inflammations may occur. These may be caused through use of the device on wounds or cuts caused by shaving, or existing wounds or ingrown hairs.
- Severe pain may occur during or after use of the device in the following cases: If you use the device on unshaven skin; if you select an intensity level that is too high for your skin tone; if you treat the same skin area several times with light pulses; if you use the device on open wounds, inflammations, tattoos, burns, etc.
- In very rare cases, permanent scarring may occur. Scarring usually occurs in the form of a flat, white lesion on the skin. Subsequent aesthetic treatment may be required to improve the appearance of the scar.

- In very rare cases, using the device can result in a blue/violet haematoma that may last 5 to 10 days. Once the haematoma fades, rusty brown discolouration of the skin (hyperpigmentation) can occur, which may be permanent.

⚠ WARNING

Warning to sunbathers before and after using the device!

Sunbathing before using the device

1. After sunbathing, wait at least 2 weeks before using the device.
2. Before using the device, check the skin for signs of sunburn. If you have sunburn that has not healed, wait until it has completely healed to use the device.
3. Then carry out a skin sensitivity test on the tanned skin to determine which intensity level is right for you.

Sunbathing after using the device

1. After using the device, wait at least 48 hours before sunbathing. Skin that is still reddened must be kept out of the sun!
2. For the first two weeks after using the device, the treated skin must be covered with clothing or sunscreen (at least SPF 30) when you are out in the sun.

Tanning with artificial light

The instructions for “Sunbathing before using the device” and “Sunbathing after using the device” also apply to tanning using artificial light and tanning lotions.

Do not use the device for 24 hours before or after using a sauna/swimming pool, as the chlorine may lead to skin irritation.

NOTICE

- Do not use the device with attachments or accessories that are not recommended by Beurer GmbH.
- During use, do not expose the device to temperatures below 5°C or above 40°C.
- Do not use the device if you see or smell smoke during use.
- Do not use the device if the fan opening on the device is cracked, is becoming detached or is missing entirely.
- After use, store the device in its original box.
- Do not expose the device to several hours of direct sunlight or UV light.

5. DEVICE DESCRIPTION

The associated drawings are shown on page 3.

- | | | | |
|---|-----------------------------|---|------------------------------|
| 1 | Light area | 6 | Trigger button |
| 2 | Skin contact sensor | 7 | Intensity display (5 levels) |
| 3 | Cooling indicator lamp | 8 | ON/OFF/intensity button |
| 4 | “Auto flash” indicator lamp | 9 | USB-C connection |
| 5 | Skin contact indicator lamp | | |

Safety features

The device also has a built-in **UV filter** that blocks harmful UV rays.

The device has a built-in **skin contact sensor** 2 to protect the eyes. It was developed so that no light pulses can be emitted when the device is held in the air. A light pulse is only emitted when the skin contact sensor 2 is in sufficient contact with the skin.

6. FIRST STEPS

6.1 Preparing for use

Before you begin using the device, you must first make appropriate preparations and then carry out a skin sensitivity test:

1. Using the hair and skin colour chart (on the back page of these instructions for use), check which intensity level is suitable for you. An “x” in the table indicates that the device is not suitable for you to use.
2. Before each application, shave the body parts that you wish to treat. Do not use a hair removal cream instead of a razor.

3. Ensure that the skin areas to be treated are shaved, dry and clean (free from any residue such as powder, cream, make-up or deodorant).

⚠ WARNING

Do not clean the skin with flammable liquids such as alcohol or acetone shortly before using the device.

- Clean the light area **1** and skin contact sensor **2** using a dry, lint-free cloth.

- Perform a skin sensitivity test on every body part that you would like to treat 48 hours before using the device fully (see chapter “6.2 Skin sensitivity test”).

6.2 Skin sensitivity test

- Connect the mains adapter to the USB-C connection **9** on the bottom of the device.
- Insert the mains adapter into a suitable socket. The ON/OFF/intensity button **8** lights up white continuously. The device is now in standby mode.
- Briefly hold down the ON/OFF/intensity button **8** to switch on the device. On the intensity display **7**, a small light (level 1) begins to illuminate continuously in white. The cooling indicator lamp **3** also begins to light up white continuously and the ventilation starts.

In the case of facial treatment, the skin sensitivity test should be performed on an area without hair, for example below the ear or on the side of the neck.

- Place the light area **1** on the skin. Ensure that the light area **1** and skin contact sensor **2** come into full contact with the skin. The skin contact indicator lamp **5** lights up as soon as there is sufficient contact with the skin.
- If the skin contact indicator lamp **5** is illuminated, press the trigger button **6** to trigger a light pulse at level 1. You will see a bright flash and hear a popping sound.
- Position the skin contact sensor **2** on another treatment area. Increase the intensity level by pressing the ON/OFF/intensity button **8** again and release a light pulse using the trigger button **6**. Do not exceed the maximum permitted/recommended intensity for your skin tone (see the hair and skin colour chart on the back page of these instructions for use).

Only increase the intensity level if the previous level is not uncomfortable for you. If intensity level 1 is uncomfortable for you, do not increase the intensity level.




- If you experience an uncomfortable feeling after the flash, reduce the intensity.
- To end the skin sensitivity test, switch off the device. To do so, press and hold the ON/OFF/intensity button **8** for 2 seconds. The intensity display **7** goes out and only the ON/OFF/intensity button **8** remains illuminated. The ventilation switches off automatically a few seconds after switching off the device. Wait 48 hours and examine the treated skin area. If the treated skin surface appears normal (no reaction or only slight reddening), you can treat the tested skin areas at the intensity level set for the skin sensitivity test.

If your treated skin surface produces a reaction (e.g. reddening or swelling), wait until the symptoms have fully subsided and then carry out another skin sensitivity test at a lower intensity. If the treated skin surface reacts again, consult your doctor.

7. USAGE

1. Connect the mains adapter to the USB-C connection **9** on the bottom of the device. Insert the mains adapter into a suitable socket.
The ON/OFF/intensity button **8** lights up continuously. The device is now in standby mode.
2. Briefly hold down the ON/OFF/intensity button **8** to switch on the device. On the intensity display **7**, a small light (level 1) begins to illuminate continuously in white. The cooling indicator lamp **3** also begins to light up white continuously and the ventilation starts.
3. If you wish to carry out an application at a higher intensity level, increase the intensity by repeatedly pressing the ON/OFF/intensity button **8** until you have reached the desired intensity.
4. Position the light area **1** on the skin and press the trigger button **6** in order to release a light pulse. Make sure that the light area **1** is in full contact with the skin and that the skin contact indicator lamp **5** is showing that there is sufficient contact. You will see a bright flash and hear a popping sound. The device will be ready for the next pulse within 1-3 seconds (depending on the intensity level).

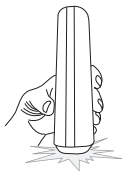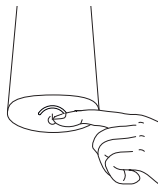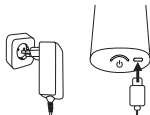7.1 Application methods

WARNING

Do not treat the same area twice within one treatment. The risk of side effects is increased if an area is treated more than once.

The pulses from the device should be applied in rows – beginning at one end of each row, continuing to the other end. This technique allows better monitoring of the areas of skin that have already been treated and helps you to avoid treating the same area more than once or treating overlapping areas of skin.

Application method 1: “Auto Flash” mode

Ensure that the skin contact sensor **2** is making sufficient contact with the skin. Press and hold the trigger button **6** for two seconds. The device automatically emits light pulses at regular intervals and the “Auto Flash” indicator lamp **4** lights up permanently. This mode does not require the trigger button **6** to be pressed manually. The “Auto Flash” mode can be ended by briefly pressing the trigger button **6**. The “Auto Flash” mode is suitable for treating large areas in particular, for example the legs.

Application method 2: Triggering individual flashes

Release the trigger button **6** after each light pulse. For the next light pulse, position the device on the area directly next to the area that has just been treated. Ensure that the light area **1** and skin contact sensor **2** come into full contact with the skin. This mode guarantees precise treatment, e.g. for knees or ankles.

If the device does not emit any light pulses and the skin contact indicator lamp **5** does not light up, this means that the skin contact sensor **2** is not making sufficient contact with the skin.

Allow the device to cool for a few minutes if it seems too hot for you. You can continue your treatment as normal after a short cooling period.

7.2 Special instructions for facial treatments

⚠ WARNING

- If the distance to the eye is less than 20cm, we recommend using safety glasses. We recommend using safety glasses approved in accordance with Regulation (EU) 2016/425 and tested in accordance with EN 166:2001 and ISO:12609-1:2021 standards.
 - Men must not use this device on the face.
 - Do not use the device on the face at an intensity level higher than 2.
 - Do not use the device on the face above the cheek line or on the eyes, eyebrows or eyelashes, as this could result in severe damage to the eyes. Do not use the device on mucous membranes in the ear or nose area.
 - The skin on your face is very sensitive and requires special care. It is essential to observe all warnings and safety notes in these instructions for use to prevent adverse consequences. We recommend that you first perform a treatment on your body and get to know the device before performing facial treatment. The device can be used to treat all areas of the face below the cheek line, such as the upper lip, cheeks, jaw, chin and neck.
1. It is essential to carry out a selective skin sensitivity test 48 hours before treatment (see chapter “6.2 Skin sensitivity test”). The skin sensitivity test should be performed on an area without hair, such as below the ear or on the side of the neck, to determine the suitable intensity level.

⚠ WARNING

Do not treat the tested area until at least 1 week after performing the skin sensitivity test.

2. Mark the treated area with a white eyeliner pencil. Use the marking as a reference line for issuing the pulses. This way you can prevent overlapping the treatment areas or treating other areas of the face by mistake.
3. Use a mirror to position the device correctly on the treatment area.
4. Follow the instructions for use and observe the following specific instructions in particular.

Treating the upper lip

Position the light area **1** of the device as shown in the figure. Trigger a pulse in every marked zone.

Avoid applying the treatment to the nostrils and lips, as these areas are more sensitive.

You can “fold in” your lips or press them together to make the treatment area more taut, or apply a thick layer of white eyeliner to the lips and around the edge of the lips. The white line ensures that the energy of the flash is reflected and does not reach your lips, should you accidentally apply a pulse to your lip.

Treating the cheeks and jaw

Position the light area **1** of the device as shown in the figure. At each point, trigger a pulse and move from one side of the cheek or jaw to the other.

Treating the neck and chin

Position the light area **1** of the device as shown in the figure. At each point, trigger a pulse and move from one side of the neck or chin to the other.

Care after facial treatment

- Apply a soothing lotion such as aloe vera.
- Do not use exfoliators, whitening lotions or similar products for at least 24 hours, as this could irritate the treated skin.
- For at least 48 hours after the treatment, avoid more than 15 minutes of direct exposure to the sun. If the treated skin is exposed to sunlight during this period, it is essential to apply sun protection with SPF 30 or higher.
- Do not remove hair on the treated areas using wax or tweezers.

Remember that it may take a complete hair growth cycle to achieve long-lasting hair removal.

7.3 Application cycle

For body treatment:

Week 1 to 2:	Use three times a week
Week 3 to 5:	Use twice a week
From week 6:	Use once or twice a month (repeat as necessary)

For facial treatment:

- The first 6 hair removal sessions using the device should be performed around two weeks apart.
- The 7th to 12th hair removal sessions using the device should be performed around four weeks apart. Afterwards, you may use the device from time to time as necessary, until long-lasting results are achieved. This corresponds to the recommended treatment plan that has proven to achieve the best results. You can, however, create a different personal treatment plan and still achieve satisfactory results.

7.4 After treatment

1. When you have completed the hair removal session, press and hold the ON/OFF/intensity button **8** for 2 seconds. The device switches to standby mode. Standby mode is indicated by the ON/OFF/intensity button **8** permanently lighting up.
2. In standby mode, the fan will continue to run for a few seconds until the device has cooled down.
3. Once the fan has switched off, remove the mains adapter from the mains socket.

⚠ CAUTION

After treatments, do not use any antiperspirants/deodorants as these may cause skin irritation.


8. CLEANING AND MAINTENANCE

We recommend cleaning your device and especially the light area **1** after each hair removal session.

When cleaning, do not use any solvents or other aggressive or abrasive cleaning agents, as the surface may become damaged.

Never immerse the device or parts of the device in water!

1. Before cleaning, remove the mains adapter of the device from the socket.
2. Clean the housing of the device with a dry, clean cloth.
3. Gently clean the light area **1** and skin contact sensor **2** using a slightly damp lint-free cloth.

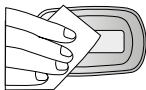9. WHAT IF THERE ARE PROBLEMS?

“The device will not switch on”

- Ensure that the mains cable is correctly connected to the device.
- Ensure that the mains adapter is inserted into a wall socket.

“No light pulse is emitted when I press the trigger button **6”**

- Ensure that the device is in full contact with the skin and that the light area **1** and skin contact sensor **2** are pressed onto the skin evenly and firmly. For your safety, the trigger button only activates a pulse if the skin contact sensor **2** has sufficient contact with the skin.
- If the skin contact indicator lamp **5** on the device is not illuminated, this indicates that the contact with the skin, as measured by the skin contact sensor **2**, is insufficient. Increase the level of contact with the skin or contact Customer Services.

10. DISPOSAL

For environmental reasons, do not dispose of the device in household waste at the end of its service life. Dispose of the device at a suitable local collection or recycling point in your country. Observe the local regulations for material disposal. Dispose of the device in accordance with EC Directive Waste Electrical and Electronic Equipment (WEEE). If you have any questions, please contact the local authorities responsible for waste disposal. You can obtain the location of collection points for old devices from the local or municipal authorities, local waste disposal companies or your retailer, for example.

11. TECHNICAL SPECIFICATIONS

Type	IPL 7800 Cool Pro	
Treatment surface	3.8 cm ²	
Technology	IPL (Intense Pulsed Light for use in the home)	
Wavelength	530 – 1100 nm	
Max. energy level	max. 5 joules/cm ²	
Intensity levels provided	5	
Time between pulses	Level 1: 1.5 s Level 2: 1.8 s Level 3: 2.1 s Level 4: 2.5 s Level 5: 2.8 s	
Maximum number of light pulses	up to 1,000,000	
Dimensions	16.4 x 6.1 x 3.9 cm	
Weight	Approx. 280 g	
Temperature	Operation	5°C to 40°C
	Storage	-20°C to 55°C
Relative humidity	Operation	≤ 80% relative humidity
	Storage	5 – 95% (non-condensing)

Mains adapter

Model number	FX48B-120400C
Input	100 – 240 V~ 50/60 Hz, 1.0 A
Output	12.0 V  4.0 A; 48 W
Protection class	This device is double-insulated (protection class II) 
Polarity	
Average efficiency during operation	$\geq 89.38\%$
No-load power consumption	$\leq 0.08 \text{ W}$

Subject to technical changes

12. WARRANTY

Further information on the warranty and warranty conditions can be found in the warranty leaflet supplied.

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Beurer GmbH • Söflinger Str. 218 • 89077 Ulm, Germany • www.beurer.com
BEURER UK Ltd., Suite 16, Stonecross Place